

MLSIGHT OPTIMIZER

POWER BOOST
WHEN NEEDED

Facilitates
near vision



 **MULTILENS**

MLSIGHT OPTIMIZER

POWER BOOST
WHEN NEEDED

Working at near distance with bifocal or progressive glasses can be strenuous, and lead to poor working position. To reduce these inconveniences and facilitate work at near distance, we have developed ML SightOptimizer.

ML SightOptimizer are spectacles with special designed optics that are easily hung on top of the habitual correction. This provides the correct individual optics, that will be better suited for near. Among the advantages are better ergonomics and more comfortable vision.

The benefits of ML SightOptimizer are an immediate “power boost”, convergence-relaxing prism, and a thin and light weight frame. ML SightOptimizer is suitable for people with reduced vision as well as for people who temporarily need extra addition.

OPTICS: Convergence optimizing lenses with a power of +1.0, +1.5, +2.0, +2.5 and +3.0

FRAME: Special designed frame, which is hung on top of habitual correction

TREATMENT: ML Prima+

WEIGHT: 15 g

FIVE FEATURES THAT ENHANCES VIEWING AT NEAR DISTANCE:

MLSIGHT OPTIMIZER

POWER BOOST

Extra addition when the occasion requires, giving better and more ergonomic seeing.

WIDE VISUAL FIELD

Open frame with clear optics throughout the lens to optimize the vision at near distance.

PRISMOPTICS

Relieving prism for a more relaxed vision at near.

ANTIREFLEX COATING

All ML SightOptimizer lenses are treated with ML Prima+. This improves the vision by creating a clearer and undisturbed image.

HANGON-FRAME

Smart and neat frame designed to fit most spectacles. The titanium frame is light and durable. The flexible temples makes it stable and comfortable for most spectacle wearers.

VARIANTS

- ML SightOptimizer 10, Increases the habitual power about 1 diopters
- ML SightOptimizer 15, Increases the habitual power about 1.5 diopters
- ML SightOptimizer 20, Increases the habitual power about 2 diopters
- ML SightOptimizer 25, Increases the habitual power about 2.5 diopters
- ML SightOptimizer 30, Increases the habitual power about 3 diopters

COST-EFFECTIVE

- The same solution works even if the habitual glasses are replaced
- Several people can use the same solution
- Lower cost than computer glasses

TEST SET 1


The test set contains following:

- 1 ML SightOptimizer 10
- 1 ML SightOptimizer 20
- 1 Fitting Guide

TEST SET 2

The test set contains following:

- 1 ML SightOptimizer 20
- 1 ML SightOptimizer 30
- 1 Fitting Guide

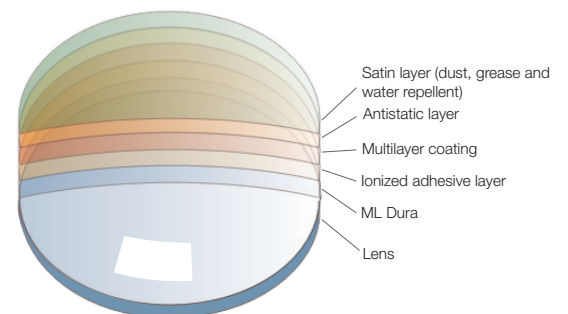

Different layers in ML Prima +



MULTILENS AB


Box 220, 435 25 Mölnlycke, Sweden

Visiting address: Skinnefjällsvägen 1

Phone: +46 (0)31 88 75 50

Fax: +46 (0)31 88 39 77

info@multilens.se, www.multilens.com



Multilens is a
unique company
within special
optics.